

Doctor's Guide to Running Your Clinic Day

Everything you need to open the system each morning, register a patient, take a measurement, and send an AI-written progress report — written for the desk, not the server room.

For the clinic front desk & consulting room

0 Before you begin

Have these five things ready before your first patient arrives. Nothing here needs to be repeated daily once it's set up.



The clinic computer

The one with SmartWeigh already installed. This is what patients' phones and the scale both talk to.



The smart scale

Eufy C1, charged and powered on, placed within about 3 metres of the computer.



Bluetooth turned on

On the clinic computer, in Windows settings — this is what lets the scale send weight readings.



The QR poster

Printed and taped up near the scale so patients can join the messaging bot in five seconds.



Internet connection

Required for AI report writing and for sending messages to patients — not required for weighing itself.

1 Starting the program

- 1 On the clinic computer's desktop, find the SmartWeigh project folder and double-click `start.bat`.
- 2 Two black windows will pop up — one titled *SmartWeigh Backend*, one titled *SmartWeigh Frontend*. Leave both open; they're doing the work in the background. You never need to type anything into them.
- 3 After a few seconds your browser opens automatically to the dashboard. If it doesn't, open Chrome or Edge and go to `localhost:3000`.

Minimize, don't close

You can minimize the two black windows once things load — just don't close them during the day. Closing either one shuts down that half of the system, and the dashboard in your browser will stop responding.

Brand-new computer only

If SmartWeigh has never been installed on this computer before, run `setup.bat` once first — it installs everything the program needs. On a computer that's already running SmartWeigh day to day, skip this; you only ever run `start.bat`.

2 Signing in

Use the clinic login on the dashboard's sign-in screen:

USERNAME	PASSWORD
doctor	changeme123

Change this before real patient data goes in

This is the default login shipped with the system. Ask your technical assistant to set a new password in the configuration file the first time SmartWeigh is installed for your clinic.

3 Registering a new patient


- 1 Click  **Patients** in the left menu.
- 2 Click **+ Register Patient** in the top-right corner.
- 3 Fill in their name, height, age, sex, and — if you've already agreed on one — their target weight.
- 4 Save. The patient now appears in your patient list and is ready for their first measurement.

Linking their Bale (messaging) account can be done now or later — see the next two sections.

4 Connecting a patient to the messaging bot

This is the step your patient does on their own phone, so you can later send them progress reports, charts, and reminders directly. It takes them under a minute and costs nothing — Bale is a free messenger app, similar to Telegram or WhatsApp.

Option A — Scan the poster (fastest)

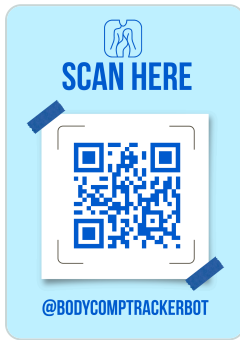

Keep this poster taped up by the scale. Ask the patient to open the **Bale** app on their phone, tap the scan icon, and point it at the code. It opens a chat with the clinic's bot immediately.

Option B — Search the bot by name

If they don't have the poster in view (or are joining from home), have them open Bale and search for:

@BodyCompTrackerBot

Type this into Bale's search bar, tap the result, then tap **Start**.

What to tell the patient


"Open the Bale app, search for BodyCompTrackerBot, and press Start. You don't need to type anything else — that just tells our system you're ready to receive your reports."

i Any message works

The patient doesn't need to type a specific word — pressing **Start** or sending any message at all is enough to register their chat with the bot. What they type isn't read or used for anything.

5 Linking their chat to the patient record

Once a patient has messaged the bot, the dashboard knows someone new has arrived — but it doesn't yet know *which* patient. You do that in one click:

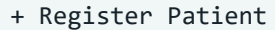- 1 Go to  Patients →  (or  Edit on an existing patient).
- 2 Look at the **Bale Account** field. If the patient has messaged the bot recently, you'll see a row of name chips: "Recent users who started the bot — click to link."
- 3 Click the chip that matches your patient's name and save. That's it — they're now linked.

✓ How you'll know it worked

A linked patient shows a small  **Bale linked** tag next to their name in the patient list, and the  Report & Send button on their profile will be able to deliver messages.

🔔 If no chip appears

Check the **Dashboard** page for a **New User** notification — it lists the chat ID Bale sent us. New-user chips only show the *last three* people who messaged the bot and haven't been linked yet, so link patients soon after they scan the code, especially in a busy clinic.

6 Taking a measurement

Open the patient's profile and go to the  Measure tab. There are three ways to record a visit — pick whichever fits the moment.

Method A — Bluetooth smart scale (recommended)

Keep the scale within 3 metres of the computer

Bluetooth Low Energy — the radio the scale uses — gets unreliable past roughly 3 metres (10 ft), and even sooner through thick walls or metal cabinets. Place the scale where the patient stands, and keep the clinic computer within that radius, or scans will time out or drop mid-reading.

- 1 Select the  Bluetooth Scale tab (it's selected by default).
- 2 Click  Start Scan. The status badge will read  Scanning... while it looks for the scale.
- 3 Ask the patient to step onto the scale **barefoot**. Bare skin against the metal foot-plates is what lets the scale measure body composition, not just weight — shoes or thick socks will give an incomplete or inaccurate reading.

- 4 Once found, the badge turns to  **Scale Connected**, then the dashboard shows a live "Waiting for patient to step on the scale..." message until the reading stabilizes.
- 5 When done, the badge reads  **Measurement Done** and weight, BMI, body fat, muscle mass, and water metrics appear on screen.
- 6 Add any clinical notes, then click  **Save Measurement**.

BADGE YOU'LL SEE	WHAT IT MEANS
 Scale Disconnected	No scan running. Click Start Scan.
 Scanning...	Looking for the scale over Bluetooth — should resolve within a few seconds if it's within range and powered on.
 Scale Connected	Found the scale — waiting for the patient's weight to stabilize.
 Measurement Done	Reading captured. Review and save.
 Bluetooth Error	Something interrupted the connection — see Troubleshooting .

Method B — Photo of another scale app (OCR)

If a patient already owns their own compatible scale (Eufy, Renpho, etc.) and shows you a screenshot from its app, you can pull the numbers straight from the photo instead of retyping them:

- 1 Select the  **OCR / Photo** tab.
- 2 Drop the screenshot in, or click to browse for it.
- 3 Click  **Extract Measurements** — the system reads the numbers off the image automatically.
- 4 **Always double-check the extracted values** before saving; photo recognition occasionally misreads a digit, especially on low-quality screenshots.
- 5 Click  **Save Measurement**.

Method C — Manual entry

For a lab-measured or verbally reported reading, select the  **Manual Entry** tab and type the values in directly.

7 Sending an AI progress report

After a measurement is saved, open the patient's profile and click

 **Report & Send** in the top-right corner.

- 1 Choose a template that fits this visit — *Comprehensive Report*, *Plateau Analysis*, *First Visit Welcome*, *Motivation Recovery*, *Lifestyle Prescription*, *Visit Reminder*, or the *6-Month Forecast Plan*. Each one asks the AI to write in a different clinical voice, in Persian, using that patient's own numbers.
- 2 Choose which extras to attach — the progress chart, the measurement table, and/or a one-page summary card.
- 3 Click  **Generate Preview** and read the AI-written message before anything is sent — you can always regenerate or edit template choice if the tone isn't right.
- 4 Click  **Send via Bale**. The patient receives it instantly in their Bale chat.

Bale account required

A patient must be linked to their Bale chat (see step 5) before anything can be sent to them. If they aren't linked yet, the button will tell you so.

8 Messaging many patients at once

For clinic-wide announcements, click  **Bulk Messaging** in the left menu.

- 1 Filter your patient list: **All Active Patients**, **Inactive 3+ Weeks**, or **Bale Linked Only**.
- 2 Pick a ready-made message — *Monthly Check-in*, *Holiday Greeting*, or *Motivational* — or write your own.
- 3 Select the patients (only those linked to Bale can actually receive it) and send.

9 Recording lab results

Open a patient's profile and select the  **Medical Tests** tab to log fasting blood sugar, HbA1c, kidney and liver panels, lipids, vitamin D, B12, and imaging notes. These are for your clinical record only — they are never sent to the patient or shown to the AI report writer.

10 Analytics

Click  **Analytics** in the left menu for a clinic-wide view — average weight loss across your patients, trend charts, and cohort comparisons. Useful for a monthly review of how the practice as a whole is doing.

11 Notifications

The  **Dashboard** keeps a running  **Recent Notifications** list with two kinds of alerts:

- **New User** — someone messaged the bot who isn't linked to a patient yet. Link them from the [Patients](#) page.
 - **Inactive** — a patient hasn't had a visit logged in a while (3 weeks by default). A good moment to send a *Visit Reminder* report.
-

12 Ending the day

Close the two black windows (*SmartWeigh Backend* and *SmartWeigh Frontend*), or simply shut the computer down — closing them stops the dashboard and the scale connection. Nothing needs to be saved manually; every measurement and message is already stored the moment you save or send it.

13 Updating the software

From time to time your developer will improve SmartWeigh — new features, fixes, refinements. Getting those updates onto the clinic computer takes two steps:

- 1 Make sure the program is **closed** (both black windows shut).
- 2 Open a terminal in the SmartWeigh project folder and run:

```
# pulls the latest version from the clinic's private repository
git pull origin master
```

- 1 Once it finishes, start the program as usual with `start.bat`.



Not confident with terminals?

This one step is the most technical part of running SmartWeigh. It's entirely reasonable to hand this single command off to whoever set the system up for you — everything else in this guide is designed to be run without any technical help at all.

If the update reports a conflict

Don't try to resolve it yourself — stop and contact your developer. This usually means a file on the clinic computer was edited directly and needs to be reconciled by hand.

14 Troubleshooting

WHAT YOU'RE SEEING	TRY THIS
Browser never opens / shows an error at localhost:3000	Check that both black windows are still open. If one crashed, close both and re-run <code>start.bat</code> .
Scale badge stuck on  Scanning... for a long time	Move the scale within 3 metres of the computer, confirm it's powered on and not asleep, and confirm Windows Bluetooth is switched on.
 Bluetooth Error appears repeatedly	Power-cycle the scale (step off, let the display go dark, step back on), then click  Start Scan again. It automatically retries in the background.
Patient can't find the bot in Bale	Confirm they typed the handle exactly: @BodyCompTrackerBot. If search fails, have them scan the poster instead.
 Send via Bale fails	Confirm the patient shows the  Bale linked tag on their profile. If linked but still failing, it's usually a temporary network issue — the system already retries automatically before giving up.
Photo upload (OCR) can't read a screenshot	Use a clearer, well-lit, non-blurry screenshot, or fall back to  Manual Entry for that visit.